



SERVICE INFORMATION LETTER

AVOX Systems Inc., 225 Erie Street, Lancaster NY 14086-9502 USA
Tel: 716-683-5100 • Fax: 716-681-1089 • Vendor Code: 53655

RECOMMENDED METHOD OF USE OF PASSENGER OXYGEN MASK (PAX MASK) AT GROUND LEVEL OR NORMAL CABIN ALTITUDE.

SUMMARY

PAX Masks manufactured by Scott Aviation or AVOX Systems are designed to provide supplemental oxygen to passengers in a loss of cabin pressure incident. They are also used to provide therapeutic or first aid oxygen in the cabins of aircraft. When used to provide supplemental oxygen at high cabin altitudes (low ambient pressure), during a loss of pressure incident, the combination of high flow of oxygen and low pressure will cause the reservoir bag of the mask to swell noticeably. However, when used to supply therapeutic or first aid oxygen at ground level or normal cabin altitudes, the bag may not noticeably inflate, even though oxygen is flowing to the mask. For the most effective use of the oxygen being supplied, the reservoir bag should be manually unfolded as directed in this Service Information Letter (SIL) when administering therapeutic or first aid oxygen.

1. PLANNING INFORMATION

A. Effectivity

(1) Aircraft

This SIL applies to all types of aircraft that have any Scott or AVOX P/N 289-601 Series, 289-701 Series, or 289-801 Series PAX Masks.

(2) Concerned equipment

This SIL applies to the PAX Masks listed in Table 1, as illustrated in Figure 1.

Table 1
End Item P/N's of Affected P/N 289 Series Passenger Oxygen Masks

289-601 Series	289-701 Series	289-801 Series
----------------	----------------	----------------

(3) Component Maintenance Manuals (CMM's)


This SIL is applicable to portions of AVOX CMM 35-22-26, CMM 35-22-27, and CMM 35-22-28. After the CMM's are revised, this SIL will not be applicable to the CMM's.

2. RECOMMENDATIONS

A. AVOX recommends that the information in this SIL be conveyed to cabin attendants who may have occasion to use PAX Masks to provide therapeutic or first aid oxygen to aircraft occupants.


(1) Prepare the oxygen supply system as instructed by your airline operating procedures.

(2) Just before or just after oxygen flow has begun, unfold PAX Mask as shown in Figure 2.

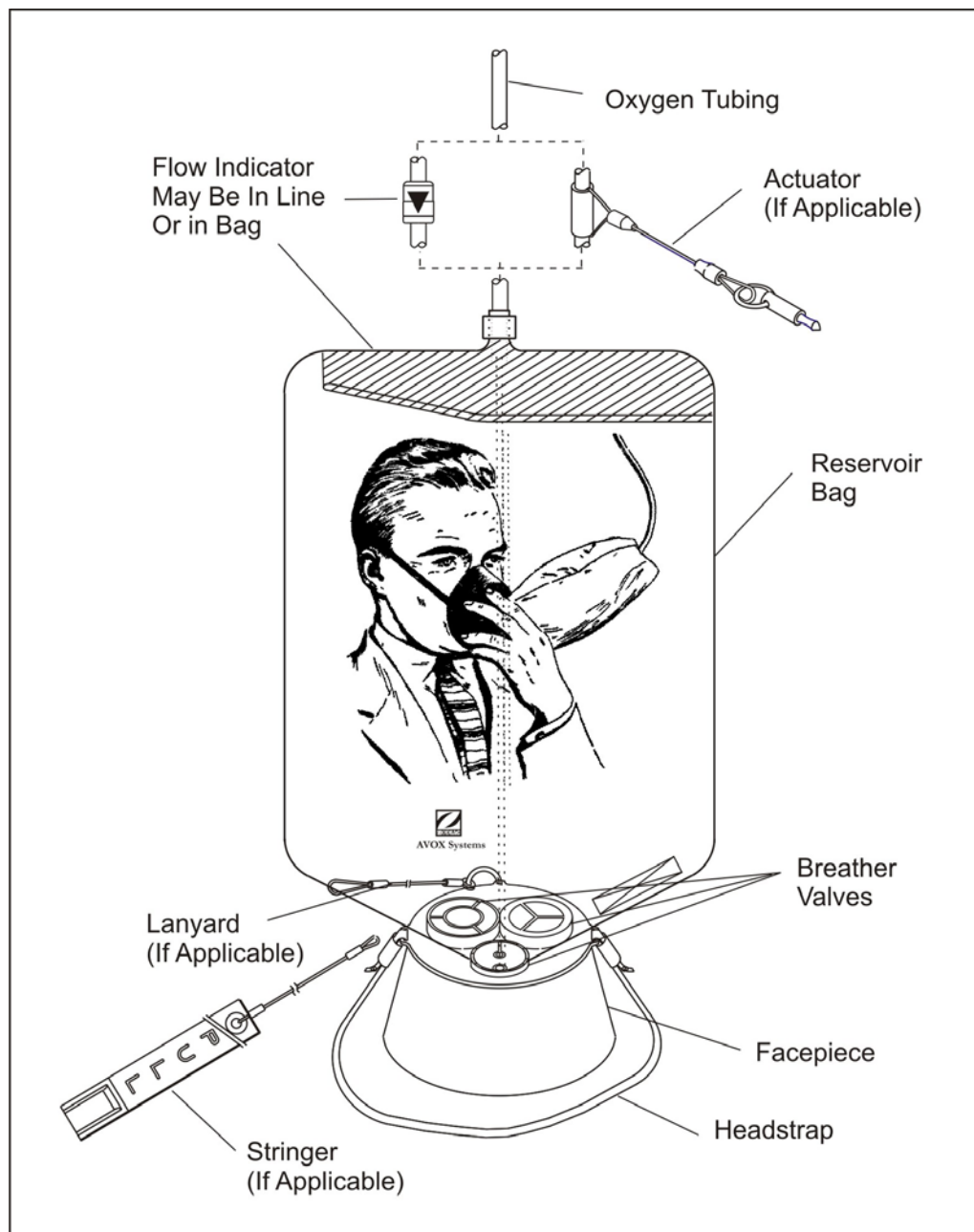

Typical PAX Mask
Figure 1

- (3) Place the mask facepiece over the passenger's mouth and nose. Use the headstrap to hold the mask in place.
- (4) Check flow indicator. Do not judge flow from the condition of the reservoir bag.

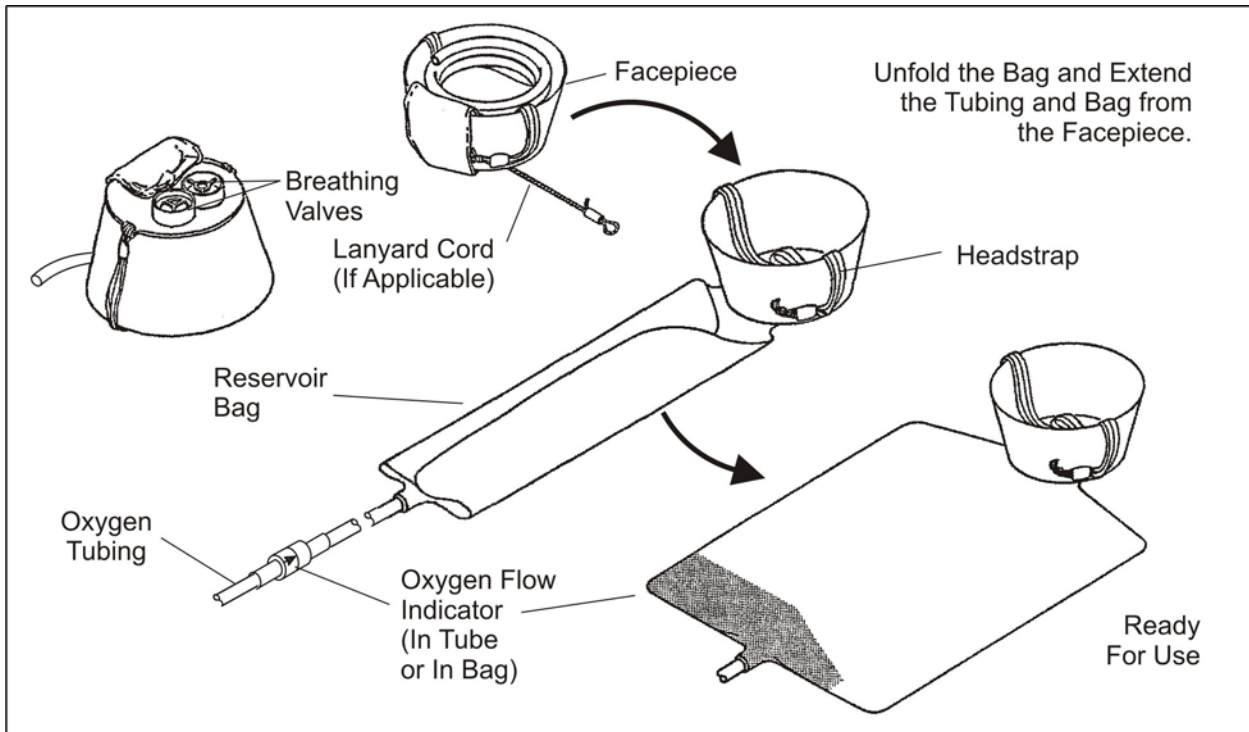

PAX Mask in Folded and Fully Unfolded Condition
Figure 2

3. APPENDIX

Contact Information:

AVOX Services Inc.

225 Erie Street
Lancaster, NY 14086-9502 USA

Tel: 716-683-5100

Fax: 716-681-1089

IN-Services Asia
80 Chun Choi Street
Tseung Kwan O Industrial Estate
New Territories
Hong Kong, China

General Manager
Sebastien WEBER
Email: weber.inasia@haeco.com
Customer Support Administration-
Shirley Huen
Email: huen.inasia@haeco.com
Phone: (852) 2260 4050
Fax: (852) 2260 4052



AVOX Systems

3. APPENDIX (Continued)

AVOX-EROS Services

6734 Valjean Ave
Van Nuys, CA 91406

Phone: 818-787-3852

Fax: 818-435-0992

Jean-Marc LeMaitre

Director of Marketing & Customer
Relations

Email: ptaylor@avox-eros.zodiac.com

818-435-0911

IN-Services Europe

7, Avenue du Parana
B.P. 303

91958 Courtaboeuf Cedex
France

Customer Services

Celine Robitaille

Phone: 33 1 64 86 69 19

Email: crobitaille@intertechnique.zodiac.com

IN-Services (Middle East)

Abu Dhabi International Airport

P.O. Box 106675

Cargo Village Area

Abu Dhabi, U.A.E.

GM and Sales Manager

Boulouafi Bouach

Phone: (971) 2 672 5257

Email: bbouach@inservices-alrumaithy.com